

4. PAEDIATRIC ANAESTHESIA EQUIPMENT

A. Endotracheal tube

UNCuffed ETT

Age	Wt (kg)	ID* (mm)	Oral length (cm)	Nasal length (cm)
Neonate	< 1	2.5	5.5	7
Neonate	1-3	3.0	6-9	7.5-11
Neonate	> 3.5	3.5	9	11
3/12	6	3.5	10	12
1 yr.	10	4.0-4.5	11	14
2 yr.	12	4.5-5.0	12	15
3 yr.	14	4.5-5.0	13	16
4 yr.	16	5.0-5.5	14	17
6 yr.	20	5.5-6.0	15	18
8 yr.	24	6.0-6.5	16	20
10 yr.	30	6.5	17	21
12 yr.	38	7.0	18	22

Although there are many guidelines available, it is the usual practice to prepare, in addition, ETTs that are half a size larger and smaller than the estimated size. The tube should not fit tightly, a small leak is advisable. Always test for leak by holding a sustained pressure of about 15 cm H₂O. After intubation, always listen in both axillae and epigastric areas.

* ID internal diameter

In general, the following formulae can be used for uncuffed ETTs in children:

ETT size (ID) : $\text{age(yrs)} / 4 + 4.5$ for patients >2 years

Insertion Depth

For children over 1 year of age:

Insertion depth (cm) for orotracheal intubation= $\text{age} / 2 + 13$

Insertion depth (cm) for nasotracheal intubation= $\text{age} / 2 + 15$

For children under 1 year of age:

Insertion depth for orotracheal intubation= $\text{weight} / 2 + 8$

Insertion depth for nasotracheal intubation= $\text{weight} / 2 + 9$

Alternatively, ETT insertion depth can be guided by the depth marker on the ETT (located at the distal end of the ETT).

The ETT should be inserted until the depth marker (dense black line) is at the level of the vocal cords. The length of the ETT at the upper incisor or nostril is then noted as the depth of insertion.

Cuffed Paediatric Endotracheal Tubes

Cuffed ETT are recommended for children >8 yrs.

Under special circumstances, we may choose to use cuffed endotracheal tubes in infants and small children. These patients often require the presence of a cuffed tube to allow adequate ventilation in the face of high airway resistance and peak ventilating pressures. It is recommended to measure cuff pressure regularly (20-30cm H₂O) We currently have cuffed ETTs from size 3-0mm ID.

In general, the following formula can be used for cuffed ETT in children:

ETT size (ID): $\text{age} / 4 + 3.5$ (except neonates)

Microcuff ETT is a specially designed paediatric tube with a high volume low pressure cuff. They afford lower risk of airway trauma and mucosal tissue injury.

Microcuff ETT Recommended Sizing for Children:

Microcuff ETT ID (mm)	Age
3.0	Term \geq 3 kg- < 8 months
3.5	8 months- < 2 years
4.0	2 years- < 4 years
4.5	4 years- < 6 years
5.0	6 years- < 8 years
5.5	8 years- < 10 years
6.0	10 years - < 12 years
6.5	12 years- < 14 years
7.0	14 years- < 16 years

B) Laryngeal Mask Airway (LMA)

Recommended sizes are:

Size	Weight (kg)	Max cuff inflation vol (ml)
1	5	up to 4
1.5	5 - 10	up to 7
2	10 - 20	up to 10
2.5	20 - 30	up to 14
3	30 - 50	up to 20
4	50 - 70	up to 30
5	70 - 100	up to 40

Cuff recommended to be inflated up to maximum of 60 cm H₂O

ProSeal LMA

ProSeal LMA is available for use in well fasted patients with no risk of regurgitation or aspiration. Choose size as for regular LMAs.

Intubation via LMA

Patients can be intubated via LMAs.

The largest ETT & FOB that can be used with each LMA is given in the table below.

LMA , FOB and ETT Size (mm)

LMA Size	Maximum FOB Size	Maximum ETT size
1	2.8mm	3.5mm Uncuffed
1.5	3.0mm	4.0mm Uncuffed
2	3.5mm	4.5mm Uncuffed
2.5	4.0mm	5.0mm Uncuffed
3	5.0mm	6.0mm Cuffed
4	5.0mm	6.0mm Cuffed
5	5.5mm	7.0mm Cuffed

*FOB : FiberOptic Bronchoscope

Our dept has FOB in the following sizes:

PAEDIATRIC ANAESTHESIA

1. 2.2mm
2. 2.8mm
3. 3.6mm

The table below indicates the smallest ETT & LMA that can be used with each FOB.

FOB	Diameter(mm))	Smallest ETT	Smallest LMA
paed	2.8	3.5*	1.0
neonatal	2.2	2.5*	1.0

* snug fit- must remove ETT blue connector before railroading over FOB. FOB must be well lubricated with silicon spray.

One Lung Ventilation in Children (Anesth analg 1999;89;1426-9)

Age(yr)	ETT(mm)	Forgarty(F)	Ardnt(F)	DLT(F)
0.5-1	3.5-4	3		
1-2	4-4.5	3		
2-4	4.5-5	3	5	
4-6	5-5.5	4-5		
6-8	5.5-6	4-5		
8-10	6 cuff	4-5	7	26
10-12	6.5 cuff	4-5		26-28
12-14	6.5-7 cuff	5-6		32

PAEDIATRIC ANAESTHESIA

14-16	7 cuff	5-6		35
16-18	7-8 cuff	7	9	35

Fogarty Catheter

Size(F)	3	4	5
Inflation vol(ml)	0.25	0.5	0.75
Diameter(mm)	8	9	10

Glidescope

Stat 0	<1.5kg
Stat 1	1.5-3.6kg
Stat2	1.8-10kg
Stat 2.5	10-28kg

C) Breathing Systems

Ayre's T- piece (Jackson Rees modification, Mapleson F)

Advantage: low resistance (no valves), minimal dead space and lightweight.

A pressure gauge should be used to measure ventilating pressures when ventilating on a T piece circuit.

Fresh gas flow (IPPV) :

1000 ml + 100 ml/kg BW/min for PaCO₂ of 4.8-5.3 kPa

or

1000 ml + 200 ml/kg BW/min

Circle system

Advantage: conserve moisture and heat, lower gas flows may be used. More economical when using expensive inhalational agents e.g. Sevoflurane. The new disposable circle system tubing and connectors are light weight.

Can be used for spontaneous respiration or IPPV. Small diameter tubing can be used for IPPV down to neonatal patients.

D) Laryngoscopes

Straight blade (Miller and Seward): use for neonates and those younger than 3 months old.

Small curved blade (Magill) can be used for those older than 3 months old.

E) Invasive Monitoring Lines

Disposable transducer sets are used – please let the Anaesthesia nurses know if you require a double or triple transducer set.

PAEDIATRIC ANAESTHESIA

The following are general recommendations based on an average sized child. If the patient is smaller or larger than expected, up or downsize as appropriate.

RADIAL Arterial lines :

Up to 3/12 or weight <5kg	:24G terumo venula
> 3/12 or weight >5kg	:22G terumo venula
> 10yr	:20G terumo venula

FEMORAL Arterial lines:

Leaderflex catheters are available in 22G 4cm for children < 20kg

Leadercath catheters are available in 20G & 18G 8 cm for children >20kg

*****Check lower limb perfusion after cannulation of femoral vessels***

Central Lines

a. Single lumen (leadercath)

Preterm neonate	<2kg	22G 4 cm
Term neonate	2-4kg	20G 8 cm
Child		20G, 18G, 8 cm


b. Double lumen

0-6 months	< 10 kg	4Fr, 5cm
------------	---------	----------

c. Triple lumen

0-6 months	<10kg	4.5Fr, 6cm
6 months - 12yrs	20-40kg	5.5Fr, 5cm, 8cm, 13cm
>12yrs	>40kg	7Fr, 16cm

Method for determining insertion depth of CVS

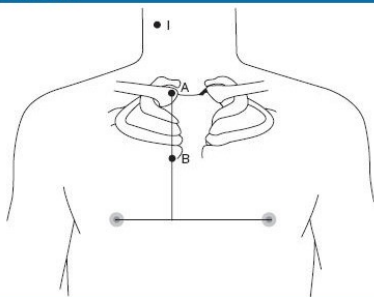

Source: Br J Anaesth © 2009 Oxford University Press

CVL tip should be positioned at junction of SVC & RA; level of carina on CXR

Two points are marked on the patient's skin during the IJV catheterization. Point A is marked at the sternal head of the right clavicle, most prominent point. Point B is marked at the midpoint of the perpendicular line from Point A to the line connecting both nipples. Point I is the insertion point of the needle. Distance from Point I to Point A and from Point A to Point B is measured. The depth of CVC is determined by adding the two measurements and subtracting 0.5 cm from this.

Practical Anatomic Landmarks for Determining the Insertion Depth of Central Venous Catheter in Paediatric Patients. H. S. Na et al. BJA 2009;102(6):820-823.

References :

- Paed Anaes. Vol 16 Issue 12 2006:1238-1243
- Black AE. Medical assessment of the paediatric patient. British Journal of Anaesthesia 1999; 83(1):3-15.
- Cook-Sather SD, Harris KA, Chiavacci R et al. A liberalised fasting guideline for formula-fed infants does not increase average gastric fluid volume before elective surgery. Anesthesia and Analgesia 2003; 96: 965-969
- Froese AB, Rose DK. A detailed analysis of T-piece systems. In: Steward, DJ (ed.) Aspects of Paediatric Anaesthesia. Amsterdam: Excerpta Medica, 1982; 101-136. Holliday MA, Segar WE. The maintenance need for water in parenteral fluid therapy. Paediatrics 1957;19:823-832 Lindhal SGE, Hulse MJ, Hatch DJ. Ventilation and gas exchange during anaesthesia and surgery in spontaneously breathing infants and children. Br J Anaesth 1984; 56: 121-129.
- Martin LD. Anesthetic Implications of an Upper Respiratory Infection in Children. Pediatric Clinics of North America 1994; 41(1) : 121-130.
- McCann ME, Kain ZN. The Management of Preoperative Anxiety in Children: An Update. Anesth Analg 2001; 93: 98-105.
- McEwan AI, Birch M, Bingham R. The preoperative management of the child with a heart murmur. Paediatric Anaesthesia 1995; 5:151-56.
- Oh TH. Formulas for calculating fluid maintenance requirements. Anesthesiology 1980; 53:351
- Practice Guidelines for Preoperative Fasting and the Use of Pharmacologic Agents to Reduce the Risk of Pulmonary Aspiration: Application to Healthy Patients Undergoing Elective Procedures. A Report by the American Society of Anesthesiologists Task Force on Preoperative Fasting. Anesthesiology 1999; 90: 896 – 905
- Rosenthal A. How to distinguish between innocent and pathologic murmurs in childhood. Pediatric Clinics of North America 1984; 31: 1229-1240.
- Tait AR, Pandit UA, Voepel-Lewis T, Munro HM, Malviya S. Use of the Laryngeal Mask Airway in Children with Upper Respiratory Tract Infections: A Comparison with Endotracheal Intubation. Anesth Analg 1998; 86: 706-711.

PAEDIATRIC ANAESTHESIA

- Tait AR, Reynolds PI, Gutstein HB. Factors that influence an anesthesiologist's decision to cancel elective surgery for the child with an upper respiratory tract infection. J Clin Anesth 1995; 7(6): 491-499.
- Van Der Walt J. Anaesthesia in children with viral respiratory tract infection. Paediatric Anaesthesia 1995; 5:287-262.
- Van der Walt JH, Moran C. An audit of perioperative management of autistic children. Paed Anaes. 2001; 11: 401-408.